Liquimmune®

LIQUID ASSAYED IMMUNOASSAY CONTROL

For In Vitro Diagnostic Use

INTENDED USE

MAS® Liquimmune® is intended for use in the clinical laboratory as a consistent test sample of known concentration for monitoring assay conditions in many immunological determinations. Include **Liquimmune®** with patient serum specimens when assaying for any of the listed constituents. Assay values are provided for the specific systems listed. The user can compare observations with expected ranges as a means of assuring consistent performance of reagent and instrument.

PRODUCT DESCRIPTION

Liquimmune® is a liquid stable control material prepared from human serum. Analyte levels are adjusted with various pure chemicals and preparations from human tissue or body fluids. Prostate specific antigen is from human seminal fluid, alpha-fetoprotein is from human cord serum and carcinoembryonic antigen is from hepatic metastases of human colon adenocarcinoma. Preservatives and stabilizers are added to maintain product integrity.

CAUTION: Liquimmune® is prepared from human source material. Components of the control which are derived from human source material have been tested using FDA accepted methods and found non-reactive for Hepatitis B Surface antigen (HBsAg), Hepatitis C (HCV), HIV-1 and HIV-2. However, no test method can offer complete assurance that products derived from human source material are free of infectious agents. This control must be handled in accordance with recommendations from Centers for Disease Control/National Institutes of Health manual, "Biosafety in Microbiological and Biomedical Laboratories, 1999."

CAUTION: The packaging of this product contains dry natural rubber.

STORAGE AND STABILITY

Unopened vials of **Liquimmune®** are stable for 90 days from receipt when stored at 2-8°C. Once opened, vials of **Liquimmune®** are stable for 30 days when stored tightly capped at 2-8°C. This product is stable until the expiration date on the box when stored at -25 to -15°C. **Self-defrosting freezers are not suitable.**

Bacterial contamination produces an increase in turbidity and/or a characteristic odor. Discard vial if evidence of microbial contamination is observed.

CONTROL RANGES

Control ranges are provided as guidelines until the laboratory has established its own statistical limits. The printed values are based upon replicate assays of representative samples by participating laboratories in accordance with established protocol and/or by direct correlation with another specific analytical system. The correlation factors used for each constituent have been derived from a comprehensive interlaboratory Quality Assurance Program. All values directly assayed have been assigned with instrument and instrument manufacturer's reagents available at the time of assay. Subsequent instrument or reagent modifications may invalidate these assigned values.

Expected values may vary slightly with different reagent and/or methodologies used. Refer to the included table for values obtained for specific systems. Values listed are specific for this lot of control only. Good laboratory practice suggests that each laboratory establish its own parameters.

INSTRUCTIONS FOR USE

Thaw control at room temperature (18-25°C) on a rocker or with periodic gentle inversion until liquid and then immediately store at 2-8°C. Thoroughly mix the contents of the vial before each use by gently inverting for several minutes. Once control is removed from 2-8°C use immediately. Open the vial and transfer the required quantity of control into a clean sample cup. Replace cap immediately and store the opened vial at 2-8°C.

ONCE THAWED, DO NOT REFREEZE THE CONTROL.

Dropper tips are available to help minimize exposure of control samples to air. DO NOT use a syringe needle to withdraw sample through the control bottle cap.

LIMITATIONS OF PROCEDURES

Compatibility of **Liquimmune®** has been demonstrated only with methods shown in this insert. Caution should be employed when using these controls with methods for which values have not been printed.

Liquimmune® contains Immunoglobulin A (IgA), Immunoglobulin G (IgG) and Immunoglobulin M (IgM); however, these constituent levels have not been adjusted but are at levels found in the human serum used in preparation of the controls. Thyroglobulin from human thyroid, Thyroid Binding Globulin from human serum, and Benzodiazepines are also weighed into **Liquimmune®**. However, no claim is made for expected values nor the stability of these constituents. Thyroglobulin and Thyroid Binding Globulin are listed here since cross reactivity with other constituents may occur with certain assay systems.

Values assigned for CK-MB were established using analytical methods measuring *mass* (ng/mL). No claim is made for expected values or stability of CK-MB determined with analytical methods measuring *activity* (IU/mL).

The specificity of antibodies used in immunoassay procedures may vary between lots. As with patient specimens, cross-reactivities between certain constituents in this control serum may occur with some reagents. Consult reagent manufacturer's package insert regarding possible interferences.

Accurate and reproducible results are dependent upon properly functioning instruments, reagents, and good laboratory technique. This product is intended for use as an assayed control for quantitative assays of listed constituents in human serum. This product is not intended for use as a calibrator.

TECHNICAL ASSISTANCE

In the USA, for technical assistance, insert updates and information, call 800-232-3342 or 510-979-5417. Outside the USA, please contact your local authorized Microgenics® distributor.

Cat. No.	Description	Size
LIG-101	Liquimmune®, Level 1	6 x 5 mL
LIG-202	Liquimmune®, Level 2	6 x 5 mL
LIG-303	Liquimmune®, Level 3	6 x 5 mL
286-606	Dropper Tips	Pkg. 100

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Architect®, AxSYM®, IMx®, TDx®, Reg.TM: Abbott Laboratories, Inc., Abbott Park, IL
 ACS-180®, Centaur, IMMUNO 1®, Reg. TM: Siemens Healthcare Diagnostics, Tarrytown, NY
 ACCESS®, BN II®, CX/LX® Series, proSPEC®, Reg. TM: Beckman Coulter, Inc. Fullerton, CA
 Mini Vidas®, Vidas®, Reg. TM: bioMérieux, Inc., Durham, NC
 Dimension® RxL, EMIT®, Stratus CS®, Reg.TM: Siemens Healthcare Diagnostics, Glasgow, DE
 Immulite®, Reg. TM: Diagnostic Products Corporation, Los Angeles, CA
 VITROS® ECI, Reg. TM: Ortho-Clinical Diagnostics, Rochester, NY
 AU Series, Reg. TM: Olympus America, Inc., Allentown, PA
 ELECSYS®, Hitachi®, Integra®, Reg.TM: Roche Diagnostics, Inc., Indianapolis, IN
 AIA®- 600, AIA®- 600-II, AIA®- Nex-IA, Reg.TM: Tosoh Medics, Inc., S. San Francisco, CA

LOT

Lot Number / Número de Lote / Chargennummer / Numéro de lot / Numero di lotto

"Use By" date/ Fecha de caducidad/ Haltbarkeitsdatum/ Date limite d'utilisation optimale/
"Utilizzare entro"Refrigerated temperature limitation / Limitación de la temperatura de refrigeración/
Kühltemperaturgrenze/ Limites de température pour un colis réfrigéré/ Limite delle temperature
di raffreddamento**REF**

Catalog Number/ Número De Catálogo/ Katalognummer/ Catalogue N°/ Numero di catalogo

Manufacturer/ Fabricante/ Hersteller/ Fabricant / Produttore

Caution/ Precaución/ Vorsicht/ Attention/ Attenzione

Consult instructions for use/ Consultar las instrucciones de uso/ Gebrauchsanweisung
beachten/ Se référer au mode d'emploi/ Consultare le istruzioni per l'uso**EC REP**Authorized Representative/ Representante autorizado/ Vertragshändler/ Mandataire/
Rappresentante autorizzato**IVD**For In Vitro Diagnostic Use/ Para uso en diagnóstico in vitro/ In-vitro-Diagnostikum/ Diagnostics
in vitro/ Per uso diagnostico in vitro**EC REP****Microgenics GmbH**

Spitalhofstrasse 94

D-94032 Passau

Germany

Tel: +49 (0) 851-88 68 90

Fax: +49 (0) 851-88 68 910

**Manufacturer / Fabricante / Hersteller / Fabricant / Produttore:****MICROGENICS CORPORATION**

46360 Fremont Blvd.

Fremont, CA 94538 USA

www.microgenics.com

ISO 13485 Certified Company

A comprehensive range of Quality Controls for your Laboratory/ Un rango comprensivo de control de calidad para su laboratorio /Eine umfassende Palette von Qualitätskontrollprodukten für Ihr Labor/ Un menu complet de contrôles de qualité pour votre laboratoire/ Range globale degli analiti per Controllo di Qualità del vostro Laboratorio

Serum Chemistry/ La Química del Suero/ Serumchemie/ Chimie de Sérum/ Chimica di Serum	
MAS[®] chemTRAK[®] · H LIQUID HUMAN BASED CHEMISTRY CONTROL/ CONTROL LÍQUIDO PARA BIOQUÍMICA/ CHEMIE-KONTROLLFLÜSSIGKEIT/ SOLUTION DE CONTROLE DE QUALITE LIQUIDE/ CONTROLLO LIQUIDO PER LA CHIMICA	MAS[®] Bilirubin LIQUID ASSAYED BILIRUBIN CONTROL/ CONTROL LÍQUIDO VALORADO PARA BILIRUBINA/ ANALYSIERTE BILIRUBIN-KONTROLLFLÜSSIGKEIT/ CONTROLE DE QUALITE LIQUIDE POUR LE DOSAGE DE BILIRUBINE/ CONTROLLO LIQUIDO TESTATO PER BILIRUBINA
MAS[®] chemTRAK[®] · Platinum LIQUID CHEMISTRY CONTROL/ CONTROL LÍQUIDO PARA BIOQUÍMICA/ CHEMIE-KONTROLLFLÜSSIGKEIT/ SOLUTION DE CONTROLE DE QUALITE LIQUIDE/ CONTROLLO LIQUIDO PER LA CHIMICA	MAS[®] Alcohol/Ammonia LIQUID ASSAYED ALCOHOL/AMMONIA CONTROL/ CONTROL LÍQUIDO VALORADO PARA ALCOHOL/AMONIACO/ ANALYSIERTE ALCOHOL/AMMONIA-KONTROLLFLÜSSIGKEIT/ CONTROLE DE QUALITE LIQUIDE POUR LE DOSAGE D'ALCOOL ET D'AMMONIAQUE/ CONTROLLO LIQUIDO TESTATO PER ALCOOL/AMMONIACA
MAS[®] Lipid LIQUID ASSAYED LIPID CONTROL/ CONTROL LÍQUIDO ANALIZADO DE LÍPIDOS/ FLUSSIGE ANALYSIERTE LIPIDKONTROLLE / CONTROLE LIQUIDE POUR DOSAGE DES LIPIDES/ CONTROLLO CAMPIONATO LIQUIDO LIPIDI	MAS[®] CSF LIQUID ASSAYED CEREBROSPINAL FLUID CONTROL/ CONTROL LÍQUIDO VALORADO PARA LÍQUIDO CEFALORRAQUÍDEO/ ANALYSIERTE LIQUOR CEREBROSPINALIS-KONTROLLFLÜSSIGKEIT/ CONTROLE DE QUALITE LIQUIDE POUR LES DOSAGES DANS LE L.C.R./ CONTROLLO LIQUIDO TESTATO PER LIQUIDO CEREBROSPINALE
Urine Chemistry/ Bioquímica de la Orina/ Urinchemie/ Dosage des Urines/ Chimica delle Urine	
MAS[®] UrichemTRAK LIQUID ASSAYED URINE CONTROL/ CONTROL LÍQUIDO VALORADO PARA ORINA/ ANALYSIERTE URIN-KONTROLLFLÜSSIGKEIT/ CONTROLE DE QUALITE LIQUIDE POUR LE DOSAGE DES URINES/ CONTROLLO LIQUIDO TESTATO PER LE URINE	MAS[®] UA Control LIQUID ASSAYED URINALYSIS CONTROL/ CONTROL LÍQUIDO VALORADO PARA ORINA/ ANALYSIERTE URIN-KONTROLLFLÜSSIGKEIT/ CONTROLE DE QUALITE LIQUIDE POUR LE DOSAGE DES URINALYSIS/ CONTROLLO LIQUIDO TESTATO PER LE URINALYSIS
MAS[®] Urinalysis Control LIQUID ASSAYED URINALYSIS CONTROL/ CONTROL LÍQUIDO VALORADO PARA ANÁLISIS DE ORINA/ ANALYSIERTE URINSTATUS-KONTROLLFLÜSSIGKEIT/ CONTROLE DE QUALITE LIQUIDE POUR LES ANALYSES D'URINES/ CONTROLLO LIQUIDO TESTATO PER ANALISI DELLE URINE	
Immunoassay/ Inmunoanálisis/ Immunoassay/ Immunologique/ Immunologici	
MAS[®] Liquimmune[®] LIQUID ASSAYED IMMUNOASSAY CONTROL/ CONTROL LÍQUIDO VALORADO PARA IMMUNOANÁLISIS/ ANALYSIERTE IMMUNOASSAY-KONTROLLFLÜSSIGKEIT/ CONTROLE DE QUALITE LIQUIDE POUR LE DOSAGE IMMUNOLOGIQUE/ CONTROLLO LIQUIDO TESTATO PER SAGGI IMMUNOLOGICI	MAS[®] T-Marker LIQUID ASSAYED TUMOR MARKER CONTROL/ CONTROL LÍQUIDO VALORADO PARA TUMOR MARKER/ ANALYSIERTE TUMOR MARKER-KONTROLLFLÜSSIGKEIT/ CONTROLE DE QUALITE LIQUIDE POUR LE MARQUEUR TUMORAUX/ CONTROLLO LIQUIDO TESTATO PER MARCATORI TUMORALI
MAS[®] PAR TDM LIQUID ASSAYED TDM CONTROL/ CONTROL LÍQUIDO VALORADO PARA FÁRMACOS/ ANALYSIERTE THERAPIEMEDIKATION-KONTROLLFLÜSSIGKEIT/ CONTROLE DE QUALITE LIQUIDE POUR LE DOSAGE DES MEDICAMENTS THERAPEUTIQUES/ CONTROLLO LIQUIDO TESTATO PER FARMACI TERAPEUTICI	
Cardiac/ Cardíacos/ Herz/ Cardiaques/ Cardiacci	
MAS[®] Cardioimmune[®] - XL LIQUID ASSAYED CARDIAC MARKER CONTROL/ CONTROL LÍQUIDO VALORADO PARA MARCADORES CARDÍACOS/ ANALYSIERTE HERZMARKER-KONTROLLFLÜSSIGKEIT/ CONTROLE DE QUALITE LIQUIDE POUR LE DOSAGE DE MARQUEURS CARDIAQUES/ CONTROLLO LIQUIDO TESTATO PER MARCATORI CARDIACI	
Protein/ Proteína/ Eiweiß/ Protéine/ Proteina	
MAS[®] Immunology LIQUID ASSAYED IMMUNOLOGY CONTROL/ CONTROL LÍQUIDO VALORADO PARA IMMUNOLOGÍA/ ANALYSIERTE IMMUNOLOGY-KONTROLLFLÜSSIGKEIT/ CONTROLE DE QUALITE LIQUIDE POUR LE DOSAGE D'IMMUNOLOGIE/ CONTROLLO LIQUIDO TESTATO PER IMMUNOLOGIA	
Specialty/ Especialidad/ Spezialität/ Spécialité/ Specialità	
MAS[®] Diabetes LIQUID ASSAYED DIABETES CONTROL/ CONTROL LÍQUIDO VALORADO PARA DIABETES / ANALYSIERTE DIABETES-KONTROLLFLÜSSIGKEIT / CONTROLE DE QUALITE LIQUIDE POUR LE DIABETE /CONTROLLO LIQUIDO TESTATO PER IL DIABETE	
Toxicology/ Toxicología/ Toxikologie/ Toxicologie/ Toxicology	
MAS[®] DOA-XSE[®] • MAS[®] DOA-XAB[®] • MAS[®] DOA GC/MS • MAS[®] Opiates 2000 • MAS[®] Opiates 2000 GC/MS LIQUID ASSAYED DRUGS OF ABUSE CONTROLS/ CONTROL LÍQUIDO VALORADO PARA DROGAS DE ABUSO/ ANALYSIERTE DROGENMISSBRAUCH-KONTROLLFLÜSSIGKEIT/ CONTROLE DE QUALITE LIQUIDE POUR LE DOSAGE DROGUES D'ABUS/ CONTROLLO LIQUIDO TESTATO PER DROGHE D'ABUSO	

Rev. 2008-0


Liquimmune®

Liquid Assayed Immunoassay Control

6 x 5mL

T.K.#59-07

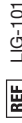LIG-101

LIA11071

2011-07-31

Microgenics Corporation
46360 Fremont Blvd.
Fremont, CA 94538 USA

CE 0123

Microgenics GmbH
Spitalhofstrasse 94
D-94032 Passau
Germany

LIG-INS-1
Rev. 31

LIG-INS
Rev. 31


Liquimmune®

Liquid Assayed Immunoassay Control

6 x 5mL

T.K.#59-09

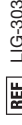LIG-303


LIA11073

2011-07-31

Microgenics Corporation
46360 Fremont Blvd.
Fremont, CA 94538 USA

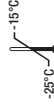CE 0123

Microgenics GmbH
Spitalhofstrasse 94
D-94032 Passau
Germany

LIG-INS-3
Rev. 31

Liquimmune®

Liquid Assayed Immunoassay Control

6 x 5mL

T.K.#59-08

LIG-202

LIA11072

2011-07-31

Microgenics Corporation
46360 Fremont Blvd.
Fremont, CA 94538 USA

CE 0123

Microgenics GmbH
Spitalhofstrasse 94
D-94032 Passau
Germany

LIG-INS-2
Rev. 31

Liquimmune®

Liquid Assayed Immunoassay Control

3 x 5mL

T.K.#59-31

LG-S1

LIA1107S

2011-07-31

Microgenics Corporation
46360 Fremont Blvd.
Fremont, CA 94538 USA

CE 0123

Microgenics GmbH
Spitalhofstrasse 94
D-94032 Passau
Germany

LIG-INS-S
Rev. 31